


RFZ-40V530

Combined Refrigerator & Freezer

SCITEK

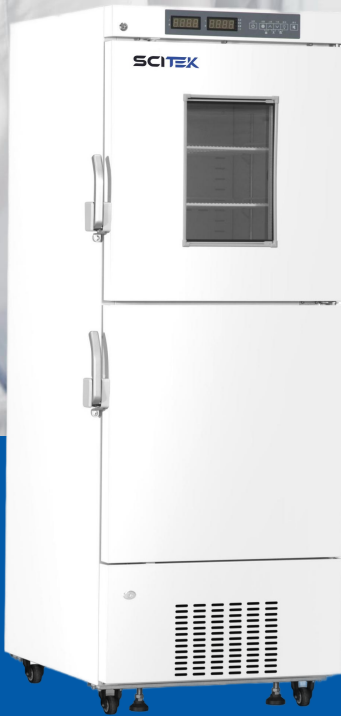

RFZ-40V530 Medical Combined Refrigerator and Freezer (upper 2~15°C and lower -40°C) is widely used in hospitals, blood stations and laboratories. The product has double chamber, double display, double control system, two sets of individual compressor and refrigeration system. Temperature in both chambers can be independently adjusted.

FEATURERS



Multi-Scene applications

Double chamber & smaller footprint



Heated glass design

Fast defogging & no condensation



Coil evaporator

Better performance, easy maintenance



Suction circulation system

Less temp. fluctuation & better



Lockable doors

Standard locks & padlocks



Large LED displays

Easy to observe and operate, built-in USB



Precise temp. control

Temp. control accurate to 0.1°C



6 types of alarm

Audible and visual alarms





Variable space

Adjustable shelves and drawers




Eco-friendly

CFC-free refrigerant and energy-saving

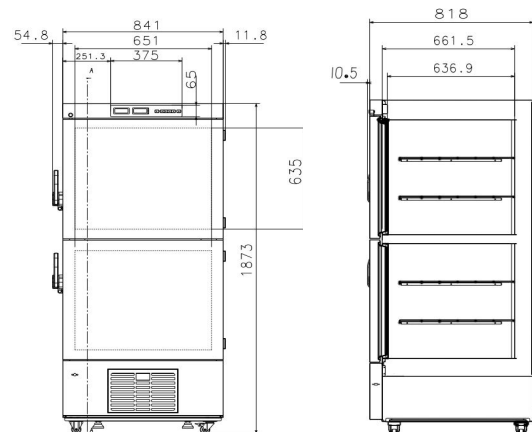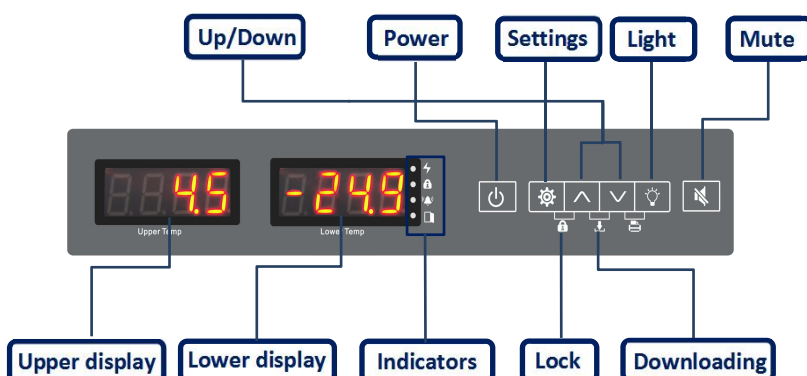

STRUCTURE

SCITEK

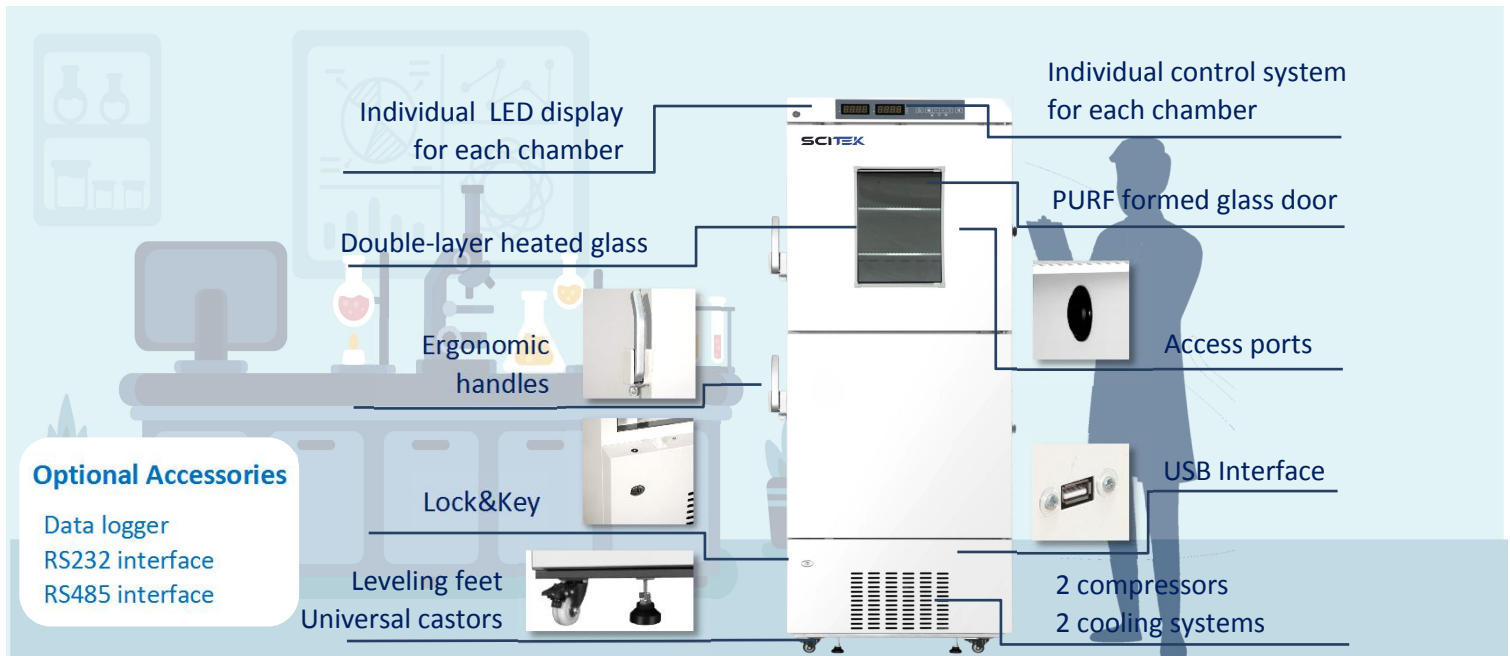

PARAMETER

Performance		Refrigeration System		Alarms	
Cooling method	Direct cooling	Temp. controller	Microprocessor	High/low temperature	Y
Defrost	Manual	Sensor	NTC	Power failure	Y
Refrigerant	R290	Display	Digital display	Sensor error	Y
Temperature range	F: $-20^\circ\text{C} \sim -40^\circ\text{C}$ R: $2^\circ\text{C} \sim 8^\circ\text{C}$	External material	SPCC sprayed steel	Controller failure	Y
Ambient temperature	$10 \sim 32^\circ\text{C}$	Condenser	Copper tube	Door ajar	Y
Noise (db)	50			Ambient temperature	Y

Electrical Data		Dimensions		Other	
Voltage frequency(V/Hz)	220V 50Hz/60Hz	Package size(W*D*H)(mm)	921×960×1950	Inside door	3 shelves & 9 drawers
Power (W)	527W	Capacity (L/cu.ft)	R:264/9.32 F:264/9.32	Remote alarm port	Optional
Current (A)	3.02A	Exterior size(W*D*H)(mm)	900×818×1870	Access port QTY/Dia.	2/25MM
Power consumption (kWh/24h)	7.06	Interior size(W*D*H)(mm)	R:765×640×635 F:651×636×635	Internal material	SPCC sprayed steel (SUS 304 for option)